



# KLT 系列激光粒度分析仪

## 平均结果报告

文件名: 钙粉 0120185.kltx

矩阵文件: Mie159\_0i00\_F100\_W.mtx

样品名: 钙粉 01

样品提供者: [REDACTED]

记录状态: Normal

记录编号: 7 到 14(8) 测量时间: 2018/5/17 20:38:09 到 2018/5/17 20:38:46

分析模式: Free Ce


分散介质: 蒸馏水

操作者: 上海盈诺 01

备注: 上海盈诺精密仪器有限公司

Dv10: 4.51 um(0.000%) Dv50: 5.79 um(0.000%) Dv90: 7.35 um(0.000%) D32: 6.01 um(0.000%)

%<45 um : 100.00% (0.000%)

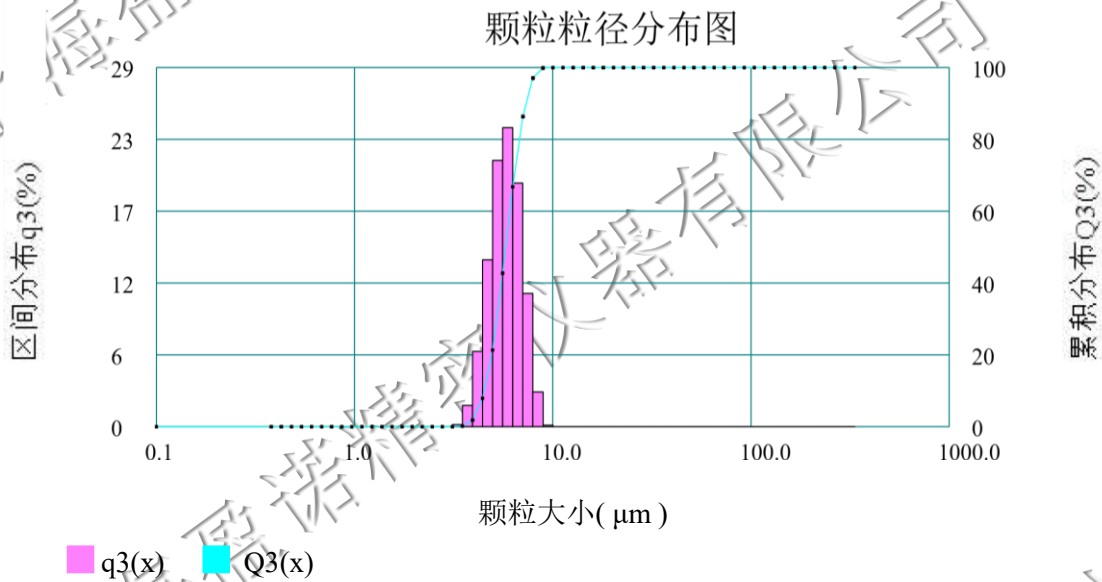

粒径(μm)	q3(x)	Q3(x)	粒径(μm)	q3(x)	Q3(x)	粒径(μm)	q3(x)	Q3(x)	粒径(μm)	q3(x)	Q3(x)
0.10	0.00	0.00	1.95	0.00	0.00	11.28	0.00	100.00	65.14	0.00	100.00
0.38	0.00	0.00	2.19	0.00	0.00	12.67	0.00	100.00	73.22	0.00	100.00
0.43	0.00	0.00	2.47	0.00	0.00	14.25	0.00	100.00	82.30	0.00	100.00
0.48	0.00	0.00	2.77	0.00	0.00	16.01	0.00	100.00	92.51	0.00	100.00
0.54	0.00	0.00	3.12	0.00	0.00	18.00	0.00	100.00	103.98	0.00	100.00
0.61	0.00	0.00	3.50	0.16	0.16	20.23	0.00	100.00	116.88	0.00	100.00
0.68	0.00	0.00	3.94	1.70	1.86	22.74	0.00	100.00	131.37	0.00	100.00
0.77	0.00	0.00	4.43	6.06	7.92	25.56	0.00	100.00	147.67	0.00	100.00
0.86	0.00	0.00	4.97	13.43	21.35	28.73	0.00	100.00	165.98	0.00	100.00
0.97	0.00	0.00	5.59	21.40	42.76	32.30	0.00	100.00	186.57	0.00	100.00
1.09	0.00	0.00	6.28	24.04	66.80	36.30	0.00	100.00	209.71	0.00	100.00
1.22	0.00	0.00	7.06	19.57	86.36	40.81	0.00	100.00	235.72	0.00	100.00
1.37	0.00	0.00	7.94	10.73	97.09	45.87	0.00	100.00	264.96	0.00	100.00
1.55	0.00	0.00	8.93	2.79	99.88	51.56	0.00	100.00	297.82	0.00	100.00
1.74	0.00	0.00	10.03	0.12	100.00	57.95	0.00	100.00	334.76	0.00	100.00